


## Perfect™ Color Protein Ladder

Cat. No.	size
E3215-01	50 loadings

**Storage Conditions:** Store at -20°C.

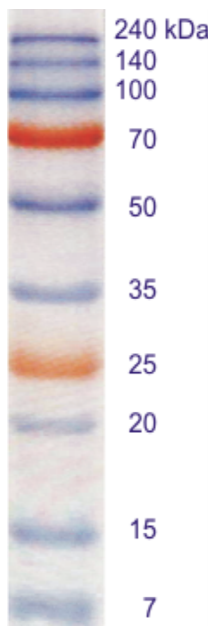

Gradient gel 8-20%

Prestained Protein Ladder, composed of precisely spaced recombinant bands, for sizing polypeptides ranging from 7 kDa to 240 kDa.

### Description:

- The Perfect™ Color Protein Ladder has ten bands corresponding to protein size: 7, 15, 20, 25, 35, 50, 70, 100, 140 and 240 kDa.
- The 25 kDa band and 70 kDa band are covalently coupled to an orange/yellow dye. The remaining bands are covalently coupled to a blue dye.
- Comes in a ready-to-use format.

### Loading:

The recommended amount of size marker to load on a mini-gel is 5 µl per lane.